


Aidal

Technical specifications

Flow	30 l/min
Vacuum	-80 kPa
Jar	1000 ml
Dimensions	35 x 14 x 28 cm
Weight	3,8 kg
Type of operation	continuous (AC)
Battery duration	40 min (depending on use)
Battery recharging time	7 h
Charging mode	electricity supply - cigarette lighter cable 12 V
Ingress Protection	IP34
Noise level at 1 meter	55,5 dBA

Contacts

 AIR LIQUIDE MEDICAL SYSTEMS S.r.l.
Via dei Prati 62
25073 Bovezzo (Brescia)
Tel. +39 030 20159.11 - Fax +39 030 2098329

www.device.airliquidehealthcare.com

Air Liquide Healthcare is a world leader in medical gases, home healthcare, hygiene products and healthcare specialty ingredients. It aims to provide customers in the continuum of care from hospital to home with medical products, specialty ingredients and services that contribute to protecting vulnerable lives.

 **Air Liquide**
HEALTHCARE

CE 0051

 **Air Liquide** | Aidal
HEALTHCARE

Homecare Suction Unit



Aidal

Homecare
Suction Unit

make it simple!

AIDAL is the brand **new suction unit** by Air Liquide Medical Systems.

The attention to each single detail makes it different: it is **studied to meet the exact needs of the user** and facilitate the care of the patient.

Every **user's interaction** with AIDAL is **therefore very simple** but even going beyond the different components, creating a full experience which is **user-friendly, effective and comfortable**.

Aidal

Totally **conceived and developed in Italy** by Air Liquide Medical Systems' team of specialists in respiratory medicine, **AIDAL is a very performing suction unit** with high flow and high vacuum pump being the best solution for any kind of users, professional or not. It is produced with **top quality materials ensuring a reliable performance and a long lasting usability.**

Care indications

- Endotracheal and bronchial aspiration
- Tracheostomy suctioning
- Nasopharyngeal and oropharyngeal aspiration



Simple to manage

- Not bulky to carry thanks to the intelligent form factors
- Comfortable and ergonomic handling

Simple to clean

- Rounded design
- Sterilizable jar included
- Single use sack compatible with the sterilizable jar

Simple to use

- Very intuitive control panel
- Easy to read vacuum gauge
- 3 Led indicators for battery status
- Vacuum knob of the ideal size to be used also with gloves
- Antibacterial filter attached to the cap of the reusable jar, improving the setting up



User friendly interface


- Easy to read vacuum gauge
- Comfortable vacuum knob for setting the vacuum



- On/Off button with light
- Green led indicator lights up when the device is connected to the mains
- 3 leds battery indicator (green, yellow, red) helps the user understand the remaining battery life and an acoustic signal warns that the battery autonomy is almost over and it is necessary to connect the device to the mains

Effective and comfortable

- The antibacterial filter is directly attached externally to the cap of the jar, ensuring a fast and simple set up of the suction unit

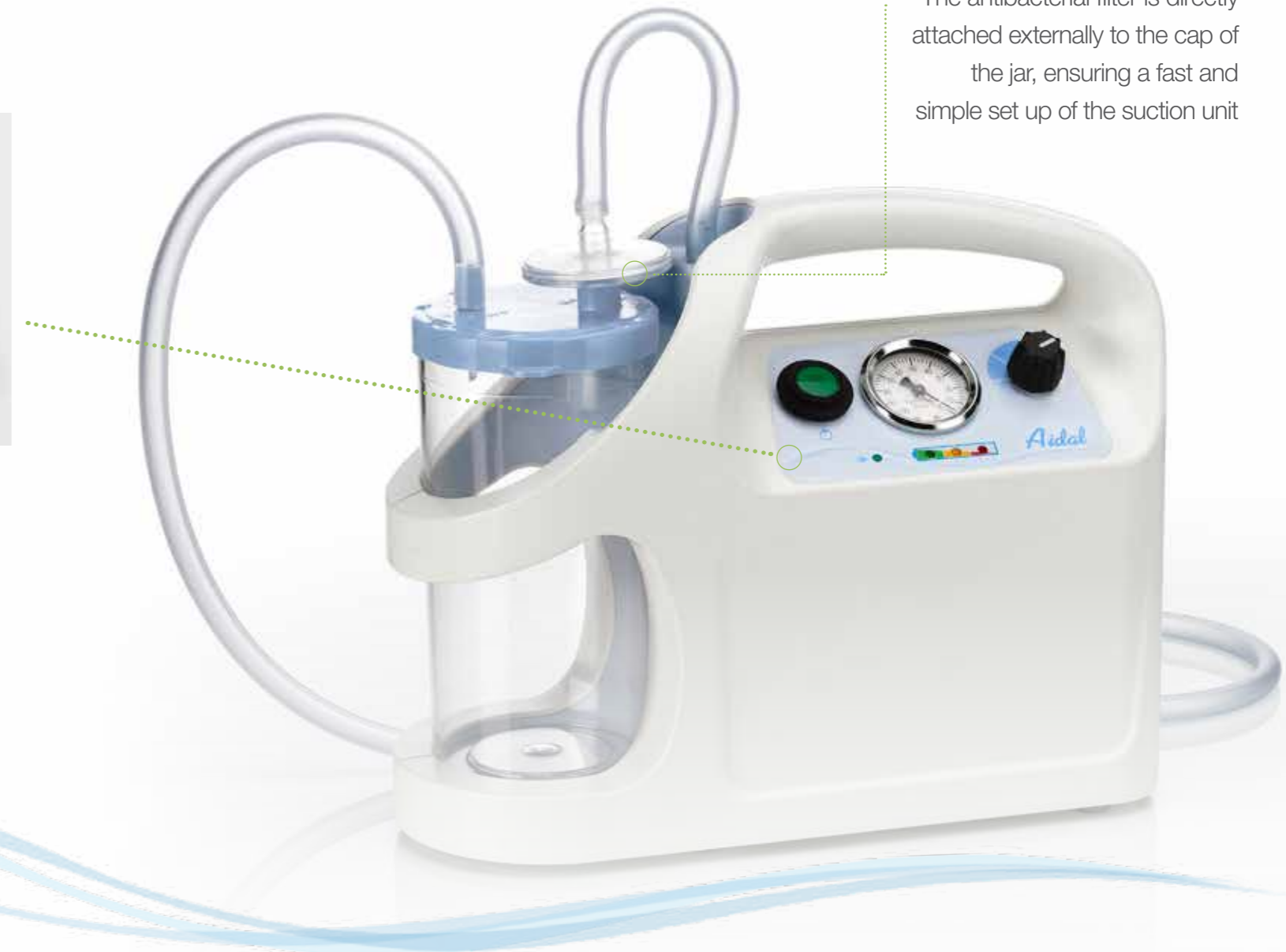

Aidal Homecare Configuration

- Suction unit
- External power supply
- Tubing
- 2 Antibacterial filters
- Reusable jar
- Fingertip controller
- Carrying bag